



NOTE: The measurements and characteristics of the products shown may be modified without notice.



BUREAU VERITAS

ECO-LABEL
ISO 14001

IBV

CE




Clinical Furniture SOFA Beds

www.medical-iberica.com

SOFA Beds

The best solution for the comfort of persons accompanying patients

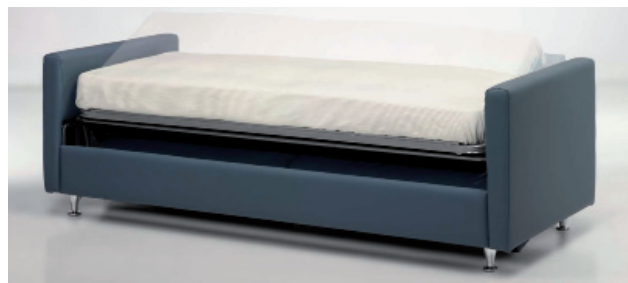

Total external measurements:	Internal measurements:
214 x 86 x 90 cm. (Height of seat 46 cm.)	193 x 57 cm. (Height of back 46 cm.)

SOFA BED MODEL MB0070

- Formed by a metal structure on a padded base.
- Divided into three sections, enabling the sofa to be put in three positions: for use as an armchair, as a chaise longue and as a horizontal bed.
- Fixed metal arms with a curved shape.
- Extensible, retractable base tray to enable it to be used as a bed or armchair.
- Back that can be converted into a bed or an armchair.
- Base formed by wooden slats measuring 80x190
- Upholstered in self-extinguishing vinyl fabric (pursuant to standard M2) which is antimicrobial, resistant to abrasion and stains, and easy to clean and maintain.

Measurements in centimetres:

Extended length: 190 cm.
Total width: 88 cm.



SOFA BED MODEL NURSE

- Sofa which can be converted into a bed using a pivoting system. Allows the mattress to be accessed without having to remove the cushions and the back.
- Extremely robust and reliable mechanism for taking out the bed, using a simple system that required no maintenance.
- Varnished steel tube structure with great resistance and durability.
- High density foam rubber mattress with great durability.
- High density foam rubber seats.
- Firm surface that prevents sinking.



SOFA BED MODEL SUMATRA

- A sofa that can be converted into a bed. Allows the mattress to be accessed by removing the cushions by means of a compass system.
- The covers can be removed, enabling the total or partial restoring of the upholstery.
- Varnished steel tube structure with great resistance and durability.
- Non-articulated welded steel mesh base that requires no maintenance.
- High density foam rubber mattress with great durability.
- The bed is always ready to use.
- High density foam rubber seats.
- Firm surface that prevents sinking.



Internal measurements:
194 cm. (Cushion interior: 64 cm.)

SOFA BED MODEL DOC

- Comprised of a wooden structure covered with polyester foam rubber.
- Interior of seat made of high density foam rubber.
- Non-deformable back with interior steel structure.
- Foldable bed with automatic mechanism for stowing cushions and back.
- Foam rubber mattress that is firm to the touch and has a padded cover.
- Detachable upholstery with Velcro fastenings.

Internal measurements:	External measurements:
114 x 48 cm. (Back: 60 cm.)	124 x 96 x 100 cm.

