

MS CEILING

Minor surgery and advanced exploration luminaire designed to be ceiling mounted. The lamp incorporates a LED light source and a highly precise articulated arm that allows the luminaire to be adjusted in any position.

The quality of the materials used in its manufacture makes it especially suitable for intensive use in surgery and examinations.

Exceptional quality of LED light

The powerful LED light source offers exceptional light quality and a lifespan of more than 50,000 hours.

The high performance of the color Ra 97% depicts the exceptional quality of white light for an excellent distinction of fabrics and its reproducing all visible colors, naturally and precisely.

Ergonomic design

The circular shape of the head and its handle in the lower part facilitate its exact positioning, fulfilling the ergonomic requirements for any use.

Designed with clean and soft shapes which facilitate hygiene and cleanliness in medical environments.

Ceiling Support


The ceiling support is included. It is manufactured specifically for the MS CEILING lamp.

The MS CEILING lamp is not compatible with any other accessory.

It is designed to be installed in ceilings up to 3 meters.

Versions

MS CEILING 45.000 lux at 50 cm
On/Off version

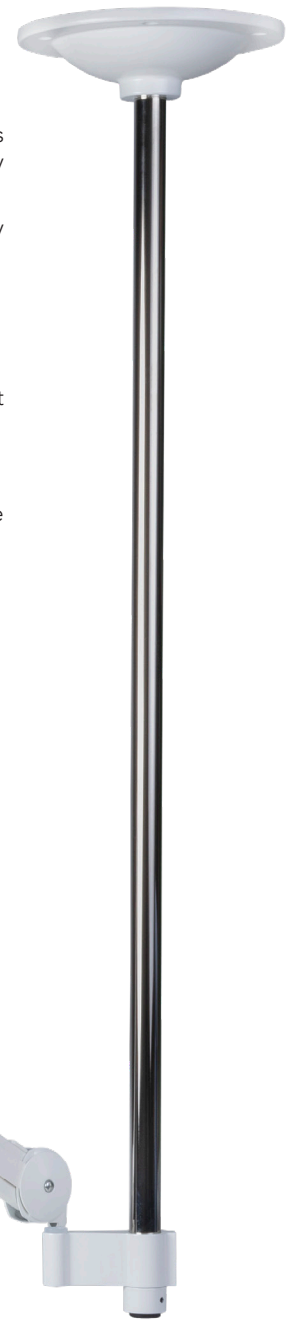

MS CEILING

Technical Details	(REF. M183500C)
Light source	LED 8W
Illumination at 50 cm	45.000 lux
Dimmable	No
Color temperature	4.500 °K
Color rendering	> 97%
Useful life	50.000 hours
Articulated Arm Length	1.050 mm
Ceiling support length	1.000 mm maximum
Color	White

