

## MT BASE: MAGNETOTHERAPY


### • 28321 GIMA MT BASE PLUS - magnetotherapy

Magneto Base Plus + is an ELF low frequency magnetotherapeutic unit designed for domiciliary applications. The most important effects of the magnetic fields are: Antiedemigen / Analgesic / Antiphlogistic / Regenerative and have applications in: traumatology, rheumatology, angiology, neurology, dermatology, endocrinology. The magnetic field produced is impulsive. The internal generator produces square waveform electrical current with frequency varying from 10 to 99 Hz and impulse time varying so to obtain a field's intensity determined from 10 to 80 Gauss. Internal sound sensor of magnetic field and safety system in case of short-circuit. Two transducers can be connected at the same time (the field produced by each transducer is half of the total).

- 28322 PLATE TRANSDUCER 70x120 mm - spare
- 28323 VELCRO ELASTIC BAND 100x8 cm - spare

#### STANDARD ACCESSORIES

Power cable	1 velcro elastic band 100xh 8 cm
2 plate transducers 80 gauss - 120x70 mm	ABS carrying case
1 velcro elastic band 60xh 8 cm	Manual (GB, IT)

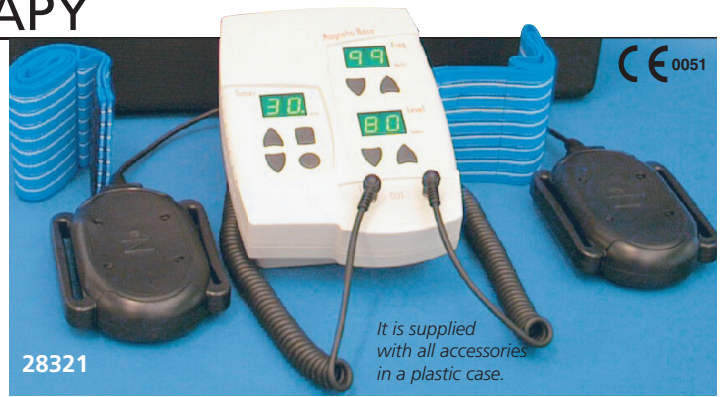

#### TECHNICAL SPECIFICATIONS

Waveform: square	Intensity of magnetic field: 80 Gauss (8mT)
Max. ratio emission/pause: 50%	Operating voltage: 230V/50-60 Hz
Selectable magnetic field from 10 to 80 Gauss	Size: 193x138xh 68 mm Weight: 1 kg
Selectable frequency: from 10 to 99 Hz	Safety class: II BF - Norms: EN60601-1 (CEI)

## IONO BASE PLUS - GALVANIC IONOPHORESIS

### • 28306 IONO BASE - PLUS - 2 channels

Local administration of medicines by ionophoretics, electric osmosis and treatment of hyperhidrosis (sweating disease) and cellulitis. Current generator for hydro galvanization, aesthetic shaping and against cellulitis slight blemish. Microprocessor controlled, with selection of parameters and internal programs. Simple and safe to use, although it maintains characteristics of a professional equipment. Made in Italy.

- 28335 RUBBER ELECTRODES 50x50 mm - box of 10
- 28336 SPONGE BAGS for electrodes 50x50 mm - box of 10

#### TECHNICAL SPECIFICATIONS

Operating Voltage: 230 V - 50/60 Hz	Size: 138x68x193 mm
Absorbed Max Power: 10 VA	Max Output Current: 20 mA
Safety class: II BF - CE 93/42/EEC. Made in Italy.	

Supplied with ABS case, velcro elastic bands, rubber conductive electrodes 80x120 mm and 50x50 mm and envelopes sponges for electrodes.




## INFRARED THERAPY LAMP - 250 W

- 28652 INFRARED THERAPY LAMP - 250 W - desk with table clamp
- 28653 INFRARED THERAPY LAMP - 250 W - trolley

The infrared light relieves the suffering caused by sore muscles and joints ankylosis. The comfortable warmth created by infrared light penetrates deep into the skin, stimulates blood circulation and warms the muscles. While the muscles are relieved, they automatically stretch and relax. At the same time, the blood circulation helps your body get rid of impurities and quickly sends oxygen-rich blood to stressed or aching muscles, causing an effective pain relief. Class II Electric.

- 28654 SPARE BULB - 250 W

Head Ø 21 cm

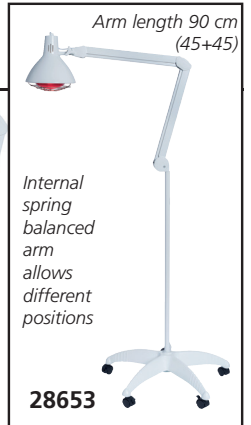

## INFRAREX: THERAPEUTIC HEATING DEVICE

- 28650 INFRAREX 230 V - 50/60 Hz
- 28651 SPARE BULB 7.5 W

Infrarex is a modern scientifically electronic therapeutic heating device which provides deep, infrared heat that can penetrate far below the surface of the skin to produce a shooting, comforting effect upon the nerves. It relieves muscular tension, reduce inflammation and provides pain relief without medication for many conditions such as: arthritis, bursitis, sinusitis, tennis elbow, back ache, neuralgia, muscular pains, aches, sprains, and other musculoskeletal aches.

**Infrarex provides pain relief without medication**, you will experience pain relief within 10-15 minutes. No concern for side effects, simple to be used.



## MOMERT INFRARED LAMP 100 W

### • 28648 MOMERT INFRARED LAMP 100 W

The lamp's deep penetrating infrared heat relaxes the muscles and stimulates the blood circulation. The pleasant heat expands the capillary veins, resulting in a better oxygen and nutrition supply to the tissues. Can be used for rheumatic diseases, allergies, cold diseases, alternative medicine (detoxification), beauty care.

Operating voltage: 230 V - 50 Hz.

- 28649 SPARE BULB - 100 W Philips

Multilanguage manual and box (GB, FR, IT, ES, PT, DE, HU, SK, LV, CZ, EE, PL, RU, BG).



Box size 185x140x225 mm

