

VISIO 2000 F.O., XENON OTO, OPHTHALMO, ENT SETS - 3.5 V

Professional diagnostic sets with fiber optic illumination system, Xenon light 3.5 V, rheostat to adjust the intensity of light. Powered by rechargeable Li-ion batteries.

All products are made of impact resistant re-inforced plastic and durable easy locking system for hospital long term use. Supplied in deluxe case, should be used with charging station 31542 and Li-ion rechargeable battery (31541). Only plug-in set (31537) as direct AC connection and need battery 31547.

- 31535 VISIO 2000 F.O. XENON OTOSCOPE - 3.5 V
- 31536 VISIO 2000 F.O. XENON OPHTHALMOSCOPE - 3.5 V
- 31537 VISIO 2000 F.O. XENON OTO-OPHTHALMOSCOPE PLUG-IN SET - 3.5 V - 2 pin plug
- 31538 VISIO 2000 F.O. XENON OTO-OPHTHALMOSCOPE - 3.5 V
- 31539 VISIO 2000 F.O. XENON DIAGNOSTIC SET - 3.5 V
- 31428 OPHTHALMOSCOPE BULB - spare
- 31544 OTOSCOPE BULB - spare
- 31541 LI-ION RECHARGEABLE BATTERY
- 31547 LI-ION RECHARGEABLE BATTERY 3.5 V PLUG-IN HANDLE for 31537

- 31542 ZOOM CHARGING STATION Supplied with EU 2 pin, UK and US plug



Swivel type window with 3X magnification for sharp images and minimal reflection



Distal F.O. illumination with super bright xenon light


Durable easy locking system

Rheostat

6 apertures to select from:

- Micro spot aperture
- Small aperture
- Large aperture
- Fixation aperture
- Slit
- Red-free filter

Enhanced Xenon super bright light

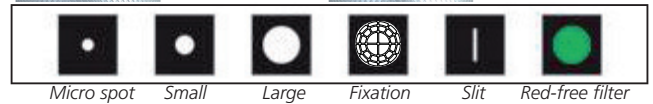

0	1	2	3	4	5	6	7	8	9	10	12	15	20	40
	1	2	3	4	5	6	7	8	9	10	15	20	25	

Colour coded 28 lenses for perfect examination (GREEN for +, RED for -)



31535

F.O. Xenon otoscope set 3.5 V
F.O. otoscope head, battery handle with bayonet lock, 1 set of 5 reusable eartips (2.5, 3.5, 4.5, 5.5 & 9.0 mm), tube of 14 disposable eartips (2.5 & 4.0 mm), spare bulb, in a rigid case.



31537

F.O. Xenon oto/ophthalmo plug-in set 3.5 V
F.O. otoscope head, F.O. ophthalmoscope head, one battery handle with bayonet lock and plug-in, 1 set of 5 reusable eartips (2.5, 3.5, 4.5, 5.5 & 9.0 mm), tube of 14 disposable eartips (2.5 & 4.0 mm), spare oto and ophthalmo bulbs, rechargeable Li-ion battery in a rigid case.



Plug-in handle



31538

F.O. Xenon oto/ophthalmo set 3.5 V
F.O. otoscope head, F.O. ophthalmoscope head, two battery handles with bayonet lock, nasal illuminator with F.O. light, 1 set of 5 reusable eartips (2.5, 3.5, 4.5, 5.5 & 9.0 mm), tube of 14 disposable eartips (2.5 & 4.0 mm), spare oto and ophthalmo bulbs in a rigid case.



31539

F.O. Xenon ENT diagnostic set 3.5 V
F.O. otoscope head, F.O. ophthalmoscope head, battery handle with bayonet lock, nasal illuminator with F.O. light, curved illuminator with bulb, plain mirror 20.0 mm, plastic tongue depressor, 1 set of 5 reusable eartips (2.5, 3.5, 4.5, 5.5 & 9.0 mm) tube of 14 disposable eartips (2.5 & 4.0 mm), spare oto and ophthalmo bulbs, in a rigid case.

GIMA VISIO 2000 WALL DIAGNOSTIC STATION - 3.5 V

- 31545 GIMA WALL DIAGNOSTIC STATION - 3.5 V
- Perfect solution for everyday clinical and hospital life. The unit has a centralized electrical supply, eliminating the need for conventional or rechargeable batteries. Wall station is equipped with two handles one with Visio 2000 otoscope Xenon light head and one with Visio 2000 ophthalmoscope Xenon light head.
- Main features:
- automatic ON/OFF switching of the handles, the instrument turns off automatically once the handle is placed in the cradle
 - non-slip ergonomic handle design
 - handle with a brightness control mechanism for light adjustment
 - spiral cable enables working at 3 metres
- Power: 115-230 V, 50-60 Hz. Disposable ear speculum see page 220



31545