


## AMBU® MARK IV RESUSCITATORS



**• 34200 AMBU® MARK IV RESUSCITATOR - adult**  
 Designed in accordance with the original double wall principle, the new 100% latex-free Ambu bags are more robust and provide users with exceptional visual and tactile lung condition feedback during manual resuscitation, making operation safer and more secure. Built-in pressure limitation prevents compression if there is any resistance in the patient's airways. Suitable for adults and children over 15 kg (3 years). Supplied with face mask N°5 and single shutter valve. Autoclavable at 134°C. Maximum volume 1300 ml.

**• 34204 AMBU® MARK IV RESUSCITATOR - baby**  
 Same as 34200, but suitable for neonates, infants and children up to 20 kg (approx 4-5 years). Supplied with face mask N°0A, oxygen reservoir tube 250 mm and pressure limiting valve 4.0 kPa (40 cm H<sub>2</sub>O). Maximum volume 300 ml.

**• 34217 AMBU MARK® IV RESUSCITATOR - adult with kit bag**  
 Same as code 34200, supplied with helicoidal mouth-wedge, tongue holding forceps, 3 Guedel airways (small, medium and large) and oxygen tube, all fitted in a resistant nylon bag.

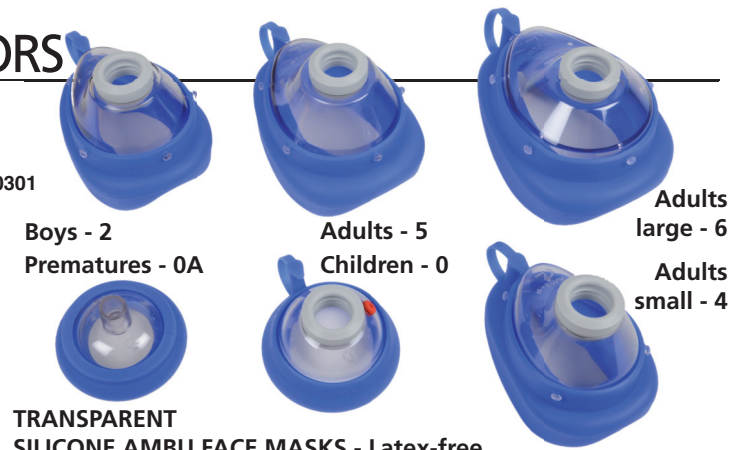

**TRANSPARENT SILICONE AMBU FACE MASKS - Latex-free**  
 The silicone face masks are self-inflating and the cuff pressure can be regulated for an easy and efficient seal. They are comfortable to grip and feature a thumb rest so rescuers can achieve a tight seal. Visual checks of bleeding, vomitus and spontaneous breathing can be made through the transparent dome. Autoclavable at 134°C

GIMA code	AMBU TRANSPARENT LATEX-FREE FACE MASKS AND RESERVOIR
34210	Ambu Baby N° 0A - newborns/prematures - blue
34205	Ambu Baby N° 0 - newborns/children - blue
34206	Ambu N° 2 - boys - blue
34207	Ambu N° 3/4 - adults/small - blue
34208	Ambu N° 5 - adults - blue
34209	Ambu N° 6 - adults/large - blue
34215	Reservoir bag 1500 ml - adult




## GIMA RESUSCITATORS - DOUBLE CHAMBER, WITH QUALITY MASK

**GIMA RESUSCITATOR BAG - with double chamber**  
 Highly professional resuscitator bag with an antistatic autoclavable and self expanding rubber internal bag and with extra soft external bag.  
 Available in 2 models: adult and infant.  
 All models are equipped with connections in accordance with EN 1281-1 directives for peep valve, intubation fitting, oxygen reservoir, oxygen intake and face masks.  
 Radiotransparent and latex free  
 All equipped with removable and autoclavable valve.  
 Made in Italy.

**• 34240 RESUSCITATOR BAG - adult**  
 Supplied with facemask N° 4 - volume 1900 cc

**• 34241 RESUSCITATOR BAG - with kit bag**  
 Same as 34240 supplied with helicoidal mouth-wedge, tongue holding forceps, 3 Guedel airways (small, medium and large) and oxygen tube, all fitted in a resistant nylon bag.

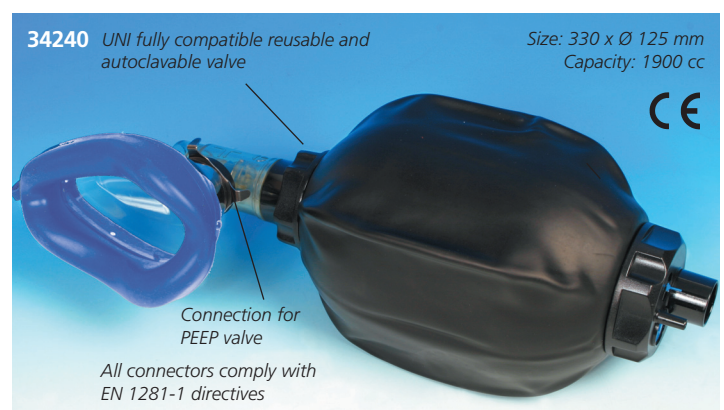

**• 34242 RESUSCITATOR BAG - infant**  
 Supplied with facemask N° 2 - volume 400 cc  
 321

