


ROLLS FOR AUTOCLAVE, STEAM AND ETO

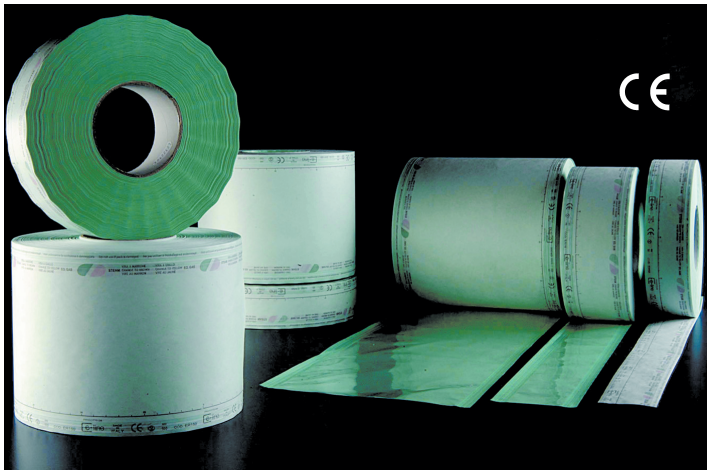

OTHER SIZES AVAILABLE

FLAT ROLLS for sterilization by steam or ethylene oxide gas
Flat rolls consist of transparent, polyester-polypropylene film coupled to medical grade paper which permit steam and gas penetration during sterilization. Steam and gas indicators provide clear colour changes when exposed to sterilization process. **Each roll is 200 m long.**

FLAT ROLLS FOR AUTOCLAVE, STEAM AND ETO				
GIMA code		External m x mm	Internal m x mm	Minimum order
35810	Flat rolls	200x50	200x36	8 rolls
35811	Flat rolls	200x75	200x53	8 rolls
35812	Flat rolls	200x100	200x78	4 rolls
35813	Flat rolls	200x150	200x128	4 rolls
35814	Flat rolls	200x200	200x178	2 rolls
35817	Flat rolls	200x250	200x228	2 rolls
35815	Flat rolls	200x300	200x278	2 rolls
35816	Flat rolls	200x400	200x378	1 roll

GUSSET ROLLS FOR AUTOCLAVE, STEAM AND ETO				
GIMA code		External m x mm x mm	Internal m x mm x mm	Minimum order
35871	Gusset rolls	100x100x25	100x78x25	4 rolls
35872	Gusset rolls	100x150x25	100x128x25	4 rolls
35873	Gusset rolls	100x200x30	100x178x30	2 rolls
35874	Gusset rolls	100x250x32.5	100x228x32.5	2 rolls
35875	Gusset rolls	100x300x40	100x278x40	1 rolls
35876	Gusset rolls	100x400x40	100x378x40	1 rolls

GUSSET ROLLS for sterilization by steam or ETO
Made with coated film PET 12/PP 40 green and 60 g medical white paper.
For sterilization of voluminous objects.
Each roll is 100 m long.

FLAT AND SELF-SEAL POUCHES FOR AUTOCLAVE



SELF-SEAL POUCHES
for sterilization by autoclave, steam or ethylene oxide gas

Self-seal pouches consist of light green, polyester-polypropylene transparent film coupled to medical grade paper which permit steam and gas penetration during sterilization.
Steam and gas indicators change colour when exposed to sterilization.

1. Reinforced corners to prevent dust from collecting near to the opening, reducing the risk of contamination during opening.
2. The laminated transparent film complies with the German standard EN 10/204. The medical paper conforms to the following standards: DIN 58953 (German), Pharmacopoeia 589 (French), pr EN 868-3 (European) and BS 6256 (English).
3. The chevrons allow rapid and correct opening.

4. These consist of one layer of laminated transparent film (special type polyester and copolymer polypropylene) and white damp-resistant medical Kraft paper, specially designed to assure a high anti-bacterial protection.

5. The heat-sealed lines (Multiline) are highly resistant and spaced equally; the external welding runs along the entire length of the pouch.

Reels and pouches are certified and are manufactured according to international standards DIN 58953 BS 6871, EN 868-5, CN 13485

OTHER SIZES AVAILABLE
minimum 10,000 pouches

FLAT POUCHES FOR AUTOCLAVE, STEAM AND ETO				
GIMA code		External mm	Internal mm	Minimum order
35880	Flat pouches	100x300	78x285	1,000
35881	Flat pouches	150x200	128x185	1,000
35882	Flat pouches	150x250	128x235	1,000
35883	Flat pouches	150x300	128x285	1,000
35884	Flat pouches	150x400	128x385	1,000
35885	Flat pouches	200x300	178x285	1,000
35886	Flat pouches	200x350	178x335	500
35887	Flat pouches	200x400	178x385	500
35888	Flat pouches	250x400	228x385	500
35889	Flat pouches	250x500	228x485	500
35890	Flat pouches	300x400	278x385	500
35891	Flat pouches	400x600	378x585	500
35895	Special size pouches			10,000

SELF-SEAL POUCHES FOR AUTOCLAVE, STEAM AND ETO				
GIMA code		External mm	Internal mm	Minimum order
35860	Self-seal pouches	90x165	68x150	1,800
35861	Self-seal pouches	90x230	68x215	1,800
35862	Self-seal pouches	140x250	118x235	1,200
35863	Self-seal pouches	190x330	168x315	1,200
35864	Self-seal pouches	190x400	168x385	1,200
35866	Self-seal pouches	300x450	278x435	400