


Personal floor scale KERN MPB

Personal floor scale with BMI function



KERN MPB-P
with stand

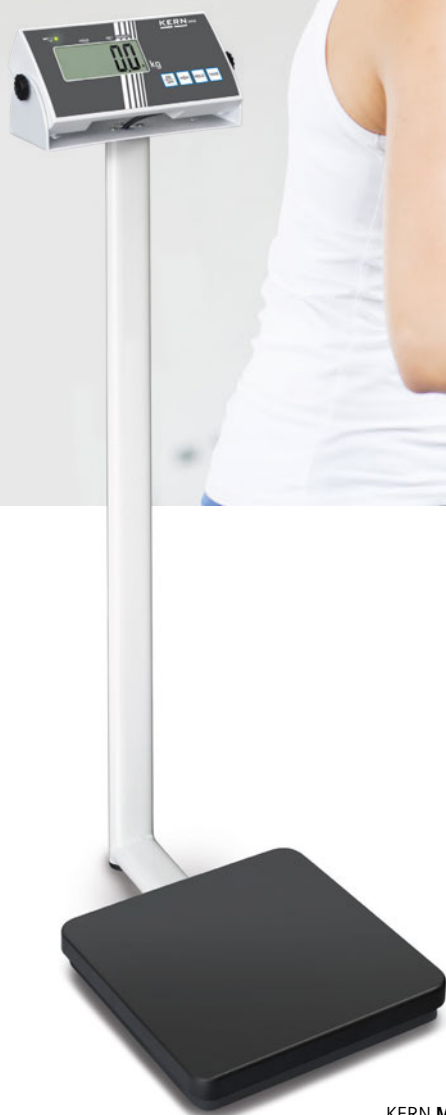

KERN MPB
with free-standing

Your KERN specialist dealer



Features

- Sturdy steel weighing plate with a non-slip and wear-resistant surface
- Easy and hygienic cleaning
- Secure and non-slip positioning with rubber feet
- Hold function:
When patients do not stand completely still, a stable weight is calculated using an average weight and this is then “frozen”. This means that you have sufficient time to attend to the patient first and then take the weight reading in peace
- Simple and convenient 4-key operation
- BMI function to determine underweight/normal weight/surplus weight
- Wall mount for the display device, only for models without a stand
- Battery- or mains-powered
- Protective working cover included with delivery

Technical data

- LCD display, digit height 25 mm
- Dimensions of weighing plate W×D×H 315×300×60 mm
- Dimensions of display device W×D×H 210×110×50 mm
- Only MPB-P: Overall dimensions (incl. stand) W×D×H 315×440×1010 mm
- Cable length of display device (only for models without a stand) approx. 120 cm
- Battery operation possible, 6×1,5 V AA, not included, operating time up to 170 h
- Mains adapter external, standard
- Net weight MPB approx. 4.6 kg, MPB-P approx. 6.4 kg

Accessories

- Mechanical Height Rod, measuring range: 60–205 cm, can be screwed onto the column from the front or the back (MPS-P) or wall-mounted, approved as a medical device according to 93/42/EEC, KERN MSF 200
- Protective working cover over the display device, scope of delivery: 5 items, KERN EOB-A02S05
- Cleaning cloths, alcohol-free cleaning cloths for wipe disinfection, quick-acting, on basis of modern quaternary ammonium compounds, and also against papova viruses. Particularly gentle on materials, and therefore very well suited for the disinfection of alcohol sensitive medical products. Fulfill the legal requirements for occupational safety in accordance with TRGS 525/540. Packaging contents 80 pcs., size 20×22 cm per cloth, KERN MYC-01

STANDARD



Model	Weighing range	Readability	
	[Max] kg	[d] kg	
KERN			
MPB 300K100	300	0,1	
MPB 300K100P	300	0,1	

Pictograms



Adjusting program CAL:

For quick setting up of the balance's accuracy. External adjusting weight required



Hold function: When patients do not stand, sit or lie completely still, a stable weight is calculated using an average weight



Battery operation:

Ready for battery operation. The battery type is specified for each device



Universal mains adapter:

with universal input and optional input socket adapters for A) EU, CH



Package shipment:

The time required for internal shipping preparations is shown in days in the pictogram